



# PICCOLO

## Lets Get Started

These instructions for use provide a guide on how to effectively use Piccolo. Please read carefully before proceeding.

### Instructions for use

1. Brush your teeth as usual using your regular toothbrush.
2. Turn on the Piccolo by pressing the base of the grip with your forefinger. You'll hear a quiet whizzing sound, which continues until you switch it off.
3. Hold the Piccolo correctly. Place the Piccolo in the palm of your hand and wrap your fingers fully and firmly around the grip. The device will not operate if held with fingertips only.
4. Adjust the tip. Squeeze and release the grip to adjust how far the tip protrudes. Release the grip to retract it.
5. Clean between the teeth. Place the tip guide gently between two teeth at the gum line, keeping it slightly away from the tooth and gum surface. This allows the tip to position itself naturally. Slowly squeeze the grip to extend the tip from the guide. The more pressure you apply, the farther the tip will move forward.
6. Retract and move on. Release the grip so the tip retracts back into the tip guide. Reposition Piccolo at the next gap.
7. Its working when. In most cases, you'll feel the tip gently massaging your gums as it moves through the spaces between your teeth. When the tip transitions from the lip side to the tongue side, you can confirm the tip is fully through the space by feeling it with your tongue against your teeth.
8. Timing matters. Insert the tip between your teeth for 1–2 seconds only—this is enough for effective cleaning. For a more thorough clean, repeat the insertion once or twice in the same space before moving to the next gap
9. For wider gaps, place the tip in the gap and gently move your hand in a small circular motion—this helps the tip contact and clean all the surrounding tooth surfaces, even in wider spaces.
10. Rinse during use. Whilst using the device, occasionally rinse the tip and tip guide under running water while the device is still on and the tip is extended. Repeat this step at the end of your cleaning session before switching off the Piccolo.
11. Need Help? [Visit the landing page](#). Review the FAQ. Contact us on 0407512378



# PICCOLO

## Lets Get Started

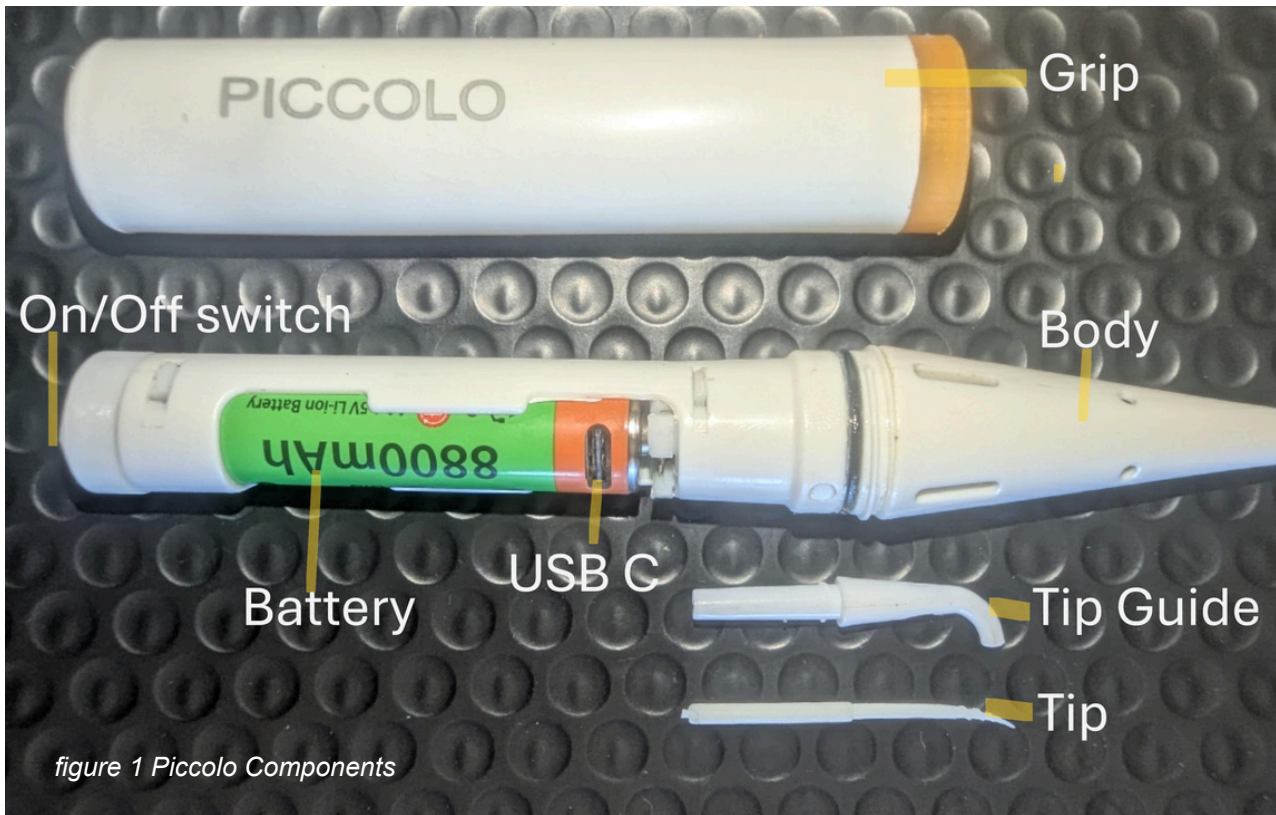


figure 1 Piccolo Components

## Warnings Contraindications & Precautions

### Warnings

- Do not force the device into tight interdental spaces. Excessive pressure with any interdental tool may cause gum irritation or disturb existing restorations.
- Do not use a damaged tip or tip guide. Replace the tip and/or tip guide if it shows signs of wear, deformation, or cracking.
- Keep out of reach of children. This device contains small components that may pose a choking hazard if detached.
- Discontinue use if you experience persistent pain, swelling, or unusual bleeding, and consult an oral health professional.

### Contraindications

Piccolo is suitable for most users. However, the following situations warrant caution or professional advice before use:

- Severe periodontal disease or active oral infections.
- Recent oral surgery or dental procedures where interdental cleaning has been temporarily restricted.
- Areas with unstable restorations, temporary crowns, loose fillings, or orthodontic appliances where mechanical vibration may not be recommended.

If you are unsure, consult your oral health professional prior to using Piccolo.



# PICCOLO

## Lets Get Started

### **Precautions**

When used as directed, Piccolo presents minimal and low-probability risks. The following precautions help ensure safe and effective use:

#### *Secure Attachment*

The Piccolo tip and tip guide are engineered to lock firmly in place. Although detachment is unlikely, always confirm that the tip guide is fully seated before use.

#### *Gentle Technique*

Piccolo is designed to clean without force.

Use light pressure and allow the device's rotating action to perform the cleaning. Avoid pushing or forcing the tip into the interdental space.

#### *Hand and Wrist Comfort*

Piccolo requires only light hand guidance.

A tight grip or unusually long sessions may cause temporary hand or wrist discomfort in sensitive individuals. Use a relaxed grip and short cleaning sessions to maintain comfort.

#### *Stop if Discomfort Occurs*

If you experience unexpected pain, sharp sensitivity, or prolonged bleeding, discontinue use and seek professional advice.

### **Potential Risks (Low Probability)**

When used according to these instructions, adverse events are rare. Potential risks include:

- Temporary gum irritation if excessive force is applied.
- Minor discomfort if the tip is used on inflamed or sensitive tissues.
- Hand or wrist fatigue with prolonged or forceful gripping.
- Tip detachment if the guide is not correctly attached (rare; preventable by checking attachment before use).